

Takiala Vave

CareSens™ NPOP

Hihitemi o te Blood Glucose Monitoring




PGJ5E3033 REV1 03/2014


1. CareSens N POP Blood Glucose Meter



2. CareSens N POP Test Strip

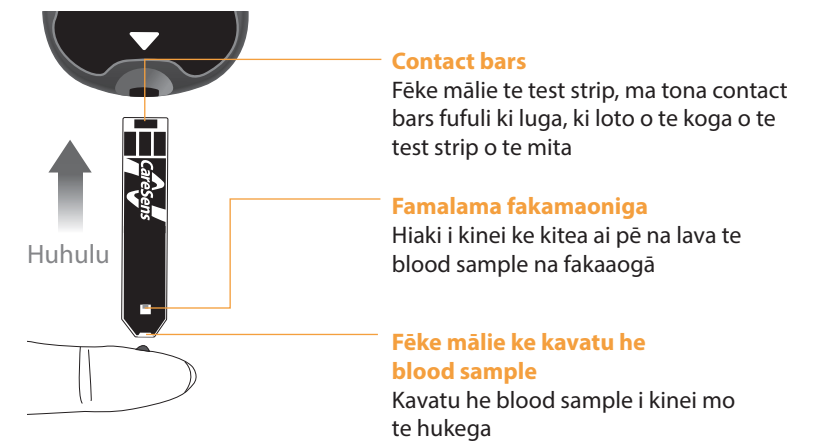

Mahini Laniheti

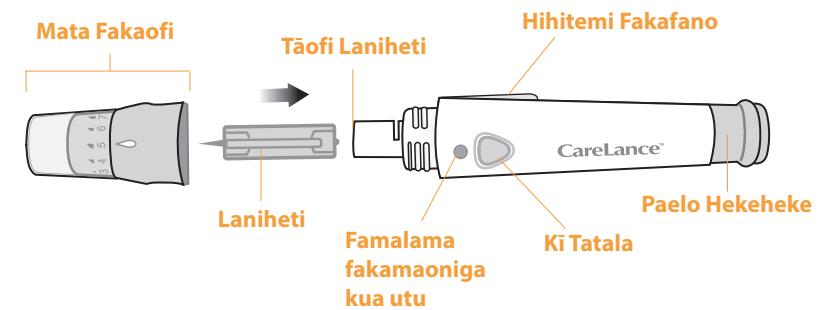

4. Hāuni te Laniheti Mahini

- Mimilo ki tua te hikulū ma kave kehe te mata fakaofi.
- Huhulu ke mau he laniheti fou ki loto o te laniheti tāofi. Kukuma mau te laniheti. Mimilo mālie ke kave kehe ai te disk puipui. Hefe te disk ke toe kāmata ai te laniheti kā uma te fakaogā. Hui te mata o te laniheti mahini.
- Filifili he loloto e manakomia i te tahi-kite-fitu (1-7) ki luga o te laniheti mahini mata fakaofi. Filifili he loloto i te fakatāvili o te vāega pito ki luga o te mata fakaofi ke pā ki te taimi kua fetau ai te numela hetu ma te fakailoga ū (arrow). Ko he hetu kāmata o te tolu(3) e fautuagia.
- Tāofi te laniheti mahini i te toho o te paelo hekeheke ke pā ki te taimi e kiliki ai te mahini. Ko te famalama e fakamaonia ai kua utu e hui ki te lanu kukula.

5. Tapenaga o te Mita ma te Test Strip

- Fufulu o lima ma te sample site. Fakamamā ma fakamago lelei.
- Huhulu he test strip ma te contact bars fufuli ki luga. Ka kitea te fakailoga ▲ i luga o te puipui (screen).

6. Kavatu te toto

- Kaumai he blood sample fakaogā ai te laniheti mahini.
- E kitea loa lava te haini ▲ i luga o te puipui (screen), kavatu loa te blood sample kite pito lauio o te test strip ke pā kite taimi e tatagi ai te mita. Kāfai ko te famalama fakamaoniga e hēki fakatumu, ka iei te fekau Er4 e kitea.
- Ko te ikuga o te huhega ka kitea kā uma te faitau e te mita āgai ki lalo mai te 5 kite 1. Ko te ikuga ka teu vave i te māfaufau o te mita.

7. Tiaki nā Laniheti kua uma te fakaogā

- Fēke ki mua te laniheti ejector i te muālima matua oi tiaki ai te laniheti kua uma te fakaogā ki he puha malu puipui lelei.

8. Fakaaliali ma Fekau

- Er-1 Ko he test strip kua uma te fakaogā na huhulu. > Toe fai te huhega i he test strip fou.
- Er-2 Ko te toto pe ko te control panel sample na kavatu muamua ka ko hēki kitea te haini ▲. > Toe fai te huhega i he test strip fou ma fakatali ke pā kite taimi e kitea ai te haini ▲ ka ko hēki kavatua te blood sample pe ko te control solution sample.
- Er-4 Ko te blood sample e mau luga kehe mai tona totō pē ko te taigole. > Toe fai te huhega kā uma te hulu o he test strip fou.
- Er-5 Ko te fekau hehē tenei e ono kitea kāfai ko te blood glucose test strip hehē e fakaogā i lō te CareSens N blood glucose test strip. > Toe fai te huhega i he CareSens N test strip.
- Er-6 Ei ei he fakafitāuli i te mita. > Nahe fakaogā te mita. Fakafehokotaki te tino fakatau fakataga o te i-SENS.



Ko te vevela i te taimi o te huhega nae i luga o te tūlaga māhani ona fakaogā. Ko te vevela i te taimi o te huhega nae i lalo o te tūlaga māhani ona fakaogā.

Kave ki he koga ei ei te vevela e māhani fakaogā (5-50°C) ma toe fai te huhega kāfai te mita ma te test strips e pā ki he vevela e māhani ona fakaogā.

Takiala Vave

