

Takiala Vave

CareSens™ N


Hihitemi o te Blood Glucose Monitoring

1. CareSens N Blood Glucose Meter

Koga o nā Fakamatalaga
Fakaaogā ke hiki ai nā fakamatalaga mai te mita ki he komipiuta kua iei he CareSens uaea

↑ Kī
Filifili pē ko te hui nā fakamatalaga

S Kī
Fakaola/tapē te mita ma ke mautinoa te filifiliga onā menu

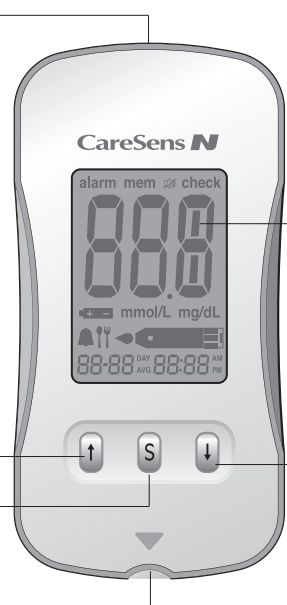


Fakaaliali
Fakailoa nā ikuga, fekau

↓ Kī
Filifili pē ko te hui nā fakamatalaga

Test Strip Port
Huhulu te test strip i kinei

Mātau: E fautuagia ke hiaki lelei te ata fakaaliali i luga o te mita e fetai tonu ma te ata fakatakitaki e vēia ona fakaali atu i luga i taimi uma e fakaola ai te mita. Nahe fakaaogā te mita pekāfai te ata fakaaliali e hē fetai tonu ma te ata fakatakitaki auā ko te mita e ono fakaali mai ni ikuga e hē hako.

2. CareSens N Test Strip

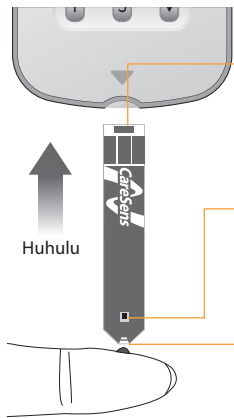


Contact bars
Fēke mālie te test strip, ma tona contact bars fufuli ki luga, ki loto o te koga o te test strip o te mita


Famalama fakamaoniga
Hiaki i kinei ke kitea ai pē na lava te blood sample na fakaaogā

Fēke mālie ke kavatu he blood sample
Kavatu he blood sample i kinei mo te huhega


3. Laniheti Mahini

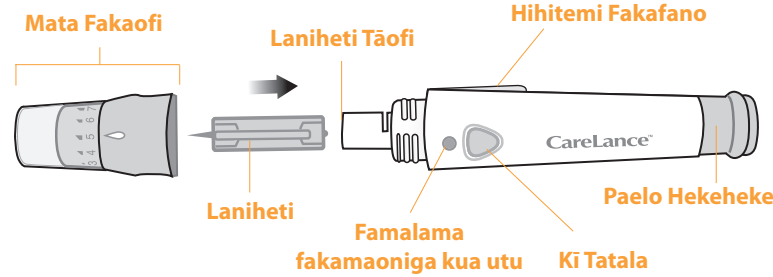


4. Hāuni te Laniheti Mahini

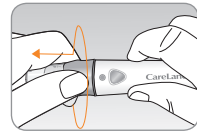


1 Mimilo ki tua te hikulū ma kave kehe te mata fakaofi.

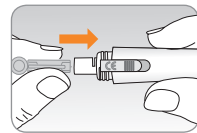

2 Huhulu ke mau he laniheti fou ki loto o te laniheti täofi. Kukuma mau te laniheti. Mimilo mālie ke kave kehe ai te disk puipui. Hefe te disk ke toe kāmata ai te laniheti kā uma te fakaaogā. Hui te mata o te laniheti mahini.

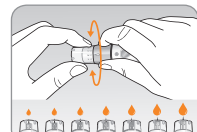

3 Filifili he loloto e manakomia i te tahi-kite-fitu (1-7) ki luga o t e laniheti mahini mata fakaofi. Filifili he loloto i te fakatāvilī o te vāega pito ki luga o te mata fakaofi ke pā ki te taimi kua fetai ai te numela heti ma te fakailoga ū (arrow). Ko he heti kāmata o te tolu(3) e fautuagia.

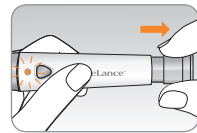

4 Täofi te laniheti mahini i te toho o te paelo hekeheke ke pā ki te taimi e kiliki ai te mahini.

5. Tapenaga o te Mita ma te Test Strip




5 Fufulu o lima ma te sample site. Fakamamā ma fakamago lelei.




6 Huhulu he test strip ma te contact bars fufuli ki luga. Ka kitea te fakailoga (screen) i luga o te ata (screen).


6. Kavatu te toto

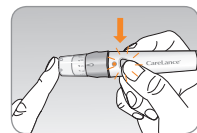

7 Kaumai he blood sample fakaaogā ai te laniheti mahini.

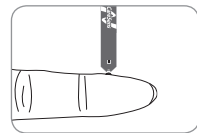


8 E kitea loa lava te haini (screen) i luga o te ata (screen), kavatu loa te blood sample kite pito lauiti o te test strip ke pā kite taimi e tatagi ai te mita. Kāfai kahē fakatumua te famalama fakamaoniga i te taimi e tatau ona fakatumua ai ona ko te kolā ma te valevalololoa (māfiafia ma te pikipikia) pe ona ko te taigole, ei ei te fekau Er4 e kitea. E fautuagia ko te faiga o te blood sample ki te test strip ke fai fakatū hako agai ki luga ki te koga e fai ai (sample site) e vēia ona fakaali mai i te ata.

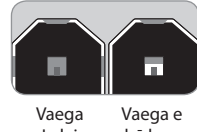


9 Ko te ikuga o te huhega ka kitea kā uma te faitau e te mita āgai ki lalo mai te 5 kite 1. Ko te ikuga ka teu vave i te māfaufau o te mita.

7. Tiaki nā Laniheti kua uma te fakaaogā

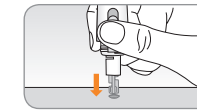

10 Tuku he tāpuni puipui ki luga o te laniheti.

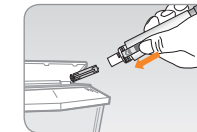

11 Fēke ki mua te laniheti ejector i te muālimamatua oi tiaki ai te laniheti kua uma te fakaaogā ki he puha malu puipui lelei.

8. Fakaaliali ma Fekau



Ko he test strip kua uma te fakaaogā na huhulu.
>Toe fai te huhega i he test strip fou.



Ko te toto pe ko te control panel sample na kavatu muamua ka ko hēki kitea te haini (screen).
>Toe fai te huhega i he test strip fou ma fakatali ke pā kite taimi e kitea ai te haini (screen) ka ko hēki kavatua te blood sample pe ko te control solution sample.



Ko te vevela i te taimi o te huhega nae i luga pe i lalo o te tūlaga māhani ona fakaaogā.
>Kave ki he koga ei ei te vevela e māhani ona fakaaogā (5-50°C/41-122°F) ma toe fai te huhega kāfai te mita ma te test strips kua pā ki he vevela e māhani ona fakaaogā.



Ko te blood sample e mau luga kehe mai tona totō pē ko te taigole.
>Toe fai te huhega kā uma te hulu o he test strip fou.




Ko te fekau hehē tenei e ono kitea kāfai ko te blood glucose test strip hehē e fakaaogā i lō te CareSens N blood glucose test strip.
>Toe fai te huhega i he CareSens N test strip.



Ei ei he fakafitāuli i te mita.
>Nahe fakaaogā te mita. Fakafehokotaki te tino fakatau fakataga o te i-SENS.



Ei ei te hehē (electronic error) e tupu i te taimi o te huhega.
>Toe fai te huhega i he test strip fou. Kāfai e kitea pea te fekau hehē (error message), fakafehokotaki tau tino fakatau fakataga o te i-SENS.

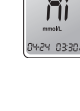

Ko te HI aikoni e kitea kāfai te lēvolo o te blood glucose e ova atu i te 33.3 mmol/L ma hypoglycemia.



Ko te Lo aikoni e kitea kāfai te blood glucose lēvolo ei lalo o te 1.1 mmol/L ma hypoglycemia.

Takiala Vave

CareSens™ **N**

